


**| We bring  
| light into the dark.**

## Instrument head for HEINE UniSpec® tubes

Illumination by fiber optic cable and projector or 3.5 V BETA handle with illumination adaptor see accessories on [page 120](#).

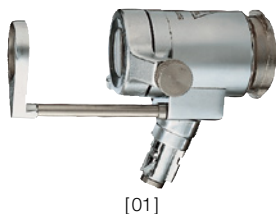

[01]

- ▶ **Annular fiber optic illumination.** Bright, reflex-free illumination
- ▶ **Magnifying lens 1.5 x.** High-resolution image
- ▶ **Hinged viewing window with insufflation port.** Improved vision by insufflation
- ▶ **Stainless steel**



[02]



[03]

### Instrument head for UniSpec disposable tubes

complete with viewing window and swivel lens 1,5 x [01]	E-003.18.100
without viewing window and swivel lens	E-003.18.098
viewing window only [02]	E-003.18.099
swivel lens 1.5 x, only (for sigmoidoscopes) [03]	E-003.18.096
same, but 1.8 x (only suitable for proctoscopes and anoscopes)	E-003.19.099

## HEINE UniSpec® disposable tubes

- ▶ **Disposable.** No cross-infection
- ▶ **No cleaning required.** Saves time
- ▶ **Moulded in patient-friendly thermoplastic.** No cold sensation on insertion
- ▶ **Light transmission through the tube wall.** Distal annular light exit. Bright, homogeneous illumination of the area examined, reflex reduction
- ▶ **Graduated inside the tube for patient comfort.** Indication of insertion length



### UniSpec disposable anoscopes (85 x 20 mm) [01]

pack of 25 disposable anoscopes	<b>E-003.19.925</b>
box of 100 (4 packs of 25) anoscopes	<b>E-003.19.911</b>

### UniSpec disposable proctoscopes (130 x 20 mm) [02]

pack of 25 disposable proctoscopes	<b>E-003.19.825</b>
box of 100 (4 packs of 25) proctoscopes	<b>E-003.19.811</b>

### UniSpec disposable sigmoidoscopes (250 x 20 mm) [03]

pack of 25 disposable sigmoidoscopes	<b>E-003.18.825</b>
box of 100 (4 packs of 25) sigmoidoscopes	<b>E-003.18.811</b>

ⓧ: Single use only.

## HEINE hygiene filter for single use

Disposable hygiene filter minimises contamination of the double bulb insufflator



disposable hygiene filters 25 pcs | pack

E-000.18.116

⊗: Single use only.

## HEINE illumination adaptor for proctological HEINE F.O. instruments

Can be attached to any HEINE 3.5 V power source to illuminate any HEINE fiber optic instrument, e.g. anoscope, proctoscope or sigmoidoscope.



Illumination adaptor

3.5V XHL

illumination adaptor

X-002.98.002

XHL xenon halogen spare bulb

X-002.88.078

► Due to local regulations, some products may not be available in your country.