



244 Bariatric Treatment Table

Safe and Efficient Care for Your Patients, You and Your Staff

850-pound Weight Capacity - With an 850-pound patient weight capacity, the Ritter 244 Barrier-Free® Bariatric Power Treatment Table is ready, willing and able to safely lift special needs

Barrier-Free® Accessibility - With an incredible 18-inch low height, the Ritter 244 is designed to safely and effectively get your special needs patients onto and off of the treatment table with minimal assistance.

The high seat height of 34 inches is designed to raise the patient to a comfortable working height, helping minimize stress and strain on your back.

The Ritter 244 features a 32-inch wide seat to help improve patient comfort and support. The seamless upholstery design looks great and offers an easy-to-clean surface.



With an ergonomic and intuitive hand control, patient positioning is simple and efficient.



The optional adjustable arm system accessory is designed to offer a sense of security for patients and a stable, secure platform to assist them in getting onto and off of the table.












18" Low Height 34" High Height



Optional casters can support 650 pounds (but should not be used to move patients).



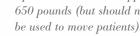



Ritter 244

S P E C I F I C A T I O N S





Fully Powered, Multi-Position Movement

Fully powered height adjustment, powered back and tilt, and angle adjustments are designed to allow for a thorough examination and safe positioning. Safety switches in the base section stop the table from moving if a slight amount of force is applied to the base cover.

Ritter 244 Bariatric Power Treatment Table

Specifications

Patient Weight Capacity: 850 lb (385 kg)

Height:

Minimum 18" (45.7 cm) Maximum 34" (86.4 cm)

Back Support Angle:

 $Minimum = 0^{\circ}$

 $Maximum = 65^{\circ}$

Seat Tilt Angle:

 $Minimum = -15^{\circ}$

 $Maximum = 15^{\circ}$ $Maximum = 15^{\circ}$

Upholstery Dimensions: Seat and Back Area: 32" W (81.3 cm)

Head and Foot Area: 30" W (76.2 cm)

Overall Length: 80" (203.2 cm)

Floor Clearance from Floor to Bottom of Base:

Leveling Feet Installed: 2.75" (7 cm)

Caster Accessory Installed: 3.5" (8.9 cm)

Paper Roll Holder Capacity: 18", 21" or 24" (45.7 cm, 53.3 cm or 61.0 cm)

Accessories

Adjustable arm system (right and left side arms included)

Casters

Foot control

3' power cord

Electrical

120 VAC, 2.5 A, 60 Hz

Safety Standards

UL 60601-1 and CAN/CSA C22.2 No. 601.1 - M90

Standard Upholstery Colors

- Mist-813
- Stone-814
- Restful Path-858
- Iron Ore-851
- Lunar Gray-845
- Obsidian-857
- Oasis-843
- Soothing Blue-855
- Dream-847
- Mineral-844
- Healing Waters-848
- Shaded Garden-853
- Citrus-849
- Sandy Retreat-850
- Curative Copper-846
- Cranberry-859
- Dark Linen-856
- Latte-860
- Robust Brown-852
- Deep Earth-854

UltraFree* Upholstery Colors

- UltraFree Stone-865
- UltraFree Restful Path-868
- UltraFree Obsidian-867
- UltraFree Harbor-861
- UltraFree Wheatgrass-862
- UltraFree Spice-863
- UltraFree Cranberry-869
- UltraFree Dark Linen-866
- UltraFree Latte-870
- UltraFree Branch-864

*Ultraleather® by Ultrafabrics

CARB 93120.2 Phase 2 Compliant

Midmark is an ISO 13485 and ISO 9001 Certified Company.

For more information or a demonstration, contact your Midmark Dealer or Midmark Corporation, Dayton, Ohio. Call: 1-800-MIDMARK Fax: 1-800-365-8631 midmark.com



© Midmark Corporation 2004 Manufactured for Midmark Corporation. Products subject to improvement changes without notice. Litho in U.S.A. 007-0537-00 Rev. H1 (1/18)