



6745 Digital Baby Scale

30 lb x 0.1 oz / 15 kg x 0.005 kg

Detecto's model 6745 Digital Baby Scale has a lock-in weight feature which compensates for typical newborn movement. When weight stabilizes, it is "locked in" and displayed. The infant can be removed and its weight is retained. The scale features a built-in battery charger, battery or AC power, a built-in measuring tape, and a large 1.0" high transfective LCD display. Serial output allows the scale to be connected to a printer or PC.



The 6745 features serial output for connectivity with a PC or printer.

ISO
REGISTERED 9001

Model 6745 Digital Baby Scale

Years of experience have gone into creating the latest in accurate, dependable infant scales... Detecto's 6745. Designed to be used with the utmost ease, this scale can weigh even the most active baby. With step-by-step instructions right on the keypad, there's never any question about what to do next. The 6745 features lock-in weight ability, built-in battery charger for NiCad or NiMH batteries (15 hours charging time), a locking pin for easily removing the tray for cleaning, and sleep and auto shutoff modes to extend battery life.

THE FLEXIBILITY YOU NEED

POWER X 3

- Rechargeable NiCad Batteries
- Standard Alkaline Batteries
- Universal AC Adapter

PRINT X 3

- P220 Thermal Label Printer
- P185 Dot Matrix Tape Printer
- Send Weight Data to a PC


6745 SPECIFICATIONS

Capacity	30 lb x 0.1 oz / 15 kg x 0.005 kg
Tray Size	26" W x 12.188" D / 66 cm W x 31 cm D
Overall Size	26" W x 14.256" D x 7.473" H / 66 cm W x 36 cm D x 19 cm H
Display	6-digit, 7-segment 1.0" / 25 mm high transfective LCD
Power	6 "C" cell batteries or universal AC adapter
Weight Units	Pounds and Ounces / Kilograms
Integrated Ruler	Measures up to 22" / 56 cm
Output	RS-232 serial
Shipping Weight	25 lb / 11 kg

6745 FEATURES

- Battery life is 250 continuous hours Alkaline / 80 continuous hours NiCad
- Sleep and Auto Shutoff modes to extend battery life
- Adjustable lock threshold
- Acceleration of gravity compensation

An integrated ruler on the baby tray measures up to 22" / 56 cm

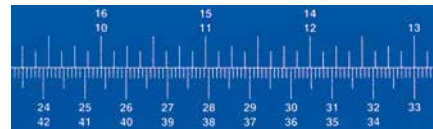

Cardinal Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.

DETECTO
A Division of Cardinal Scale Manufacturing Co.
203 E. Daugherty, Webb City, MO 64870 USA
Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-5001
www.Detecto.com

Sold By: