



Portable Digital **FLOOR SCALES**

DR400C / DR550C



BATTERY OR AC POWER

Detecto
Model DR400C
with Low-Profile
Platform and Rubber Mat

- 2.2-inch/5.6-cm high, low-profile platform for easy portability
- Coiled cable from scale base allows versatile positioning of indicator
- Slip-resistant rubber mat for traction and patient stability
- 550 lb x 0.2 lb and 400 lb x 0.5 lb models available

- Lightweight, mobile weighing platform for easy transport
- Auto Shut-Off after 2 minutes conserves battery power
- Simple four-key operation is easy to use even for novices
- Pushbutton Tare will zero out the load on the platform
- Pounds and kilograms conversion with Units key
- RS232 serial port on DR400C for EMR/EHR-ready output

DETECTO's model DR400C and DR550C digital portable floor scales will provide you with years of accurate and reliable weighing from the brand you trust in clinical weighing. The scales arrive assembled for use.

The portable lightweight platform and large bold readouts make this scale ideal for mobile clinics, home healthcare nurses, doctor's offices, general purpose use, and home bathroom weighing.

The scale's moveable display lets you place it on a work surface or mounted to the wall, wherever it's most convenient for you to view. An optional black carrying case is available.



DETECTO is the largest medical scale manufacturer in the United States providing clinical-grade accuracy and dependability for over a century.



Detecto
Model DR550C
with Stainless Steel
Low-Profile Platform



The economical DR400C and DR550C are battery powered, so you can take the scales anywhere you need them. The bold LCD is viewable in almost any lighting condition.

Versatile Display Options

The indicator's versatile mounting bracket and coiled cable allows the display to be positioned on the floor or tabletop for quick mobile weighing use or you may permanently mount it on the wall using the included screws (see below).



FLOOR / TABLETOP



WALL MOUNT

MOBILE Positioning




The bracket can temporarily rest on any flat surface for mobile weighing

FIXED Positioning



2 screws and 2 drywall anchors are included for the mounting bracket

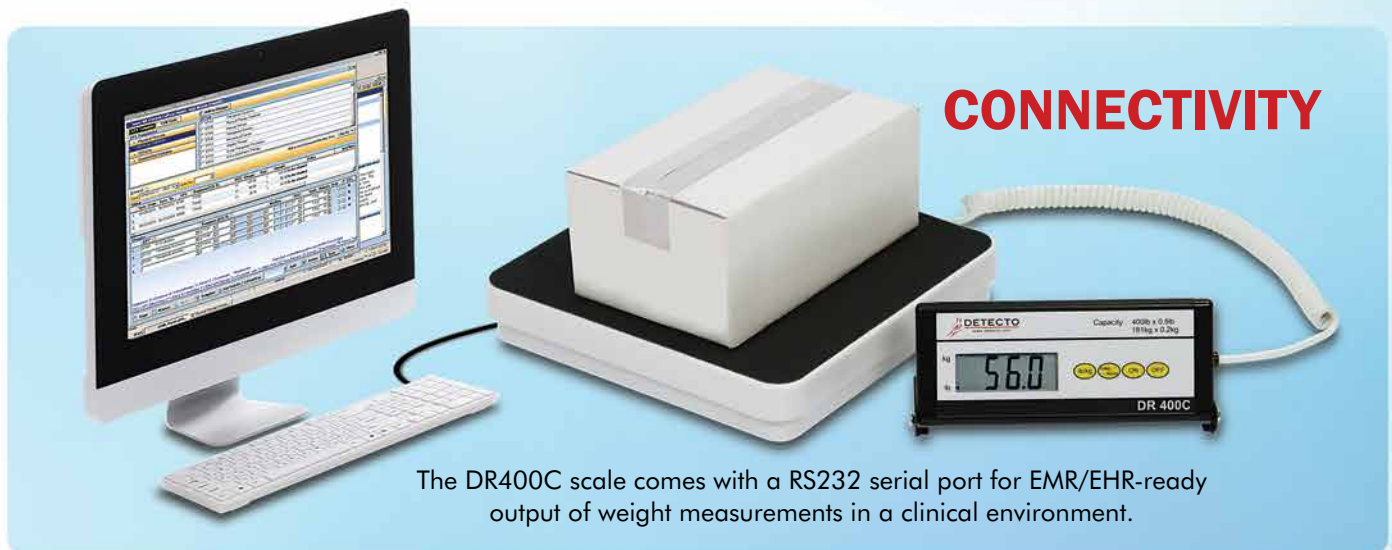

CONNECTIVITY

The DR400C scale comes with a RS232 serial port for EMR/EHR-ready output of weight measurements in a clinical environment.

Battery or AC Power

The DR400C and DR550C come with a UL listed 9V 100mA DC adapter (included) or may be powered by 6 AA batteries (not included) for mobile use.



The DR400C features a RS232 serial port (not included on DR550C).



Easy access to 6 AA batteries located conveniently underneath the platform cover.

Portable Digital Floor Scales



	DR400C	DR550C
Capacity	400 lb x 0.5 lb / 180 kg x 0.2 kg	550 lb x 0.2 lb / 250 kg x 0.1 kg
Scale Base Platform Dimensions	11.9 in W x 12.3 in D x 2.2 in H / 30.2 cm W x 31.3 cm D x 5.6 cm H	
Weight Indicator Dimensions	8.8 in W x 3.6 in D x 1.8 in H / 22.2 cm W x 9.2 cm D x 4.4 cm H (minus bracket)	
Scale Base Platform Construction	Mild steel with rubberized mat	Stainless steel with rubberized mat
Display Type	Liquid Crystal Display, 4 digit, 7 segment	
LCD Display Dimensions	2.9 in W x 1.1 in H / 7.5 cm W x 2.7 cm H overall LCD size	
LCD Digit Dimensions	0.9 in H x 0.4 in W / 2.4 cm H x 1.0 cm W per digit	
Cable Length (indicator to base)	Coiled: 33.4 in / 84.8 cm Fully Extended: 117.5 in / 298.5 cm	
Keypad	LB/KG, TARE/HOLD, ON, OFF	ON/OFF, NET/GROSS, UNIT, TARE
Connectivity Port	RS232 serial	
Power	6 AA 1.5V batteries (not included) or UL listed 9V 100mA DC adapter (included)	
Tare	100% of full scale capacity	
Zero	Established on power up routine and maintained by auto-zero circuitry	
Operating Temperature	40° F - 105° F (5° C - 40° C)	
Operating Humidity	25% ~ 95% RH	
UPC Code	809161143904	809161144000
Net Weight	10 lb / 4.5 kg	
Shipping Weight	13 lb / 5.9 kg	

Detecto Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.

OPTIONAL ACCESSORIES



Wall-Mount Grab Bar

Stainless steel wall-mount support handles (ordered 1 bar per model)


Model GBSS18-WM



Carrying Case

Black carrying case with handle and zipper for DR400C and DR550C

Model DR-CASE

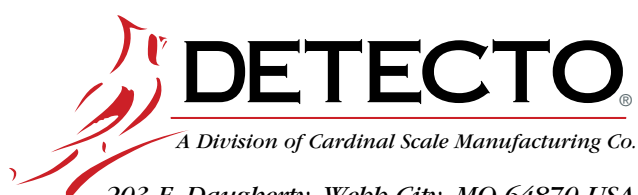

A Division of Cardinal Scale Manufacturing Co.

203 E. Daugherty, Webb City, MO 64870 USA

Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-5001

www.Detecto.com

Mobile: m.detecto.com

SOLD BY: