


**It's sure to increase
your heartbeat.
Out of pure joy.**

HEINE sphygmomanometers and stethoscopes

SHOCKPROOF

- Shockproof double-walled construction with unique manometer protection. [01]

The robust double-walled construction with a unique triple manometer mount gives ideal protection against shocks. Meets the requirements of european standard EN ISO 81060-1 and ANSI | AAMI | ISO 81060-1

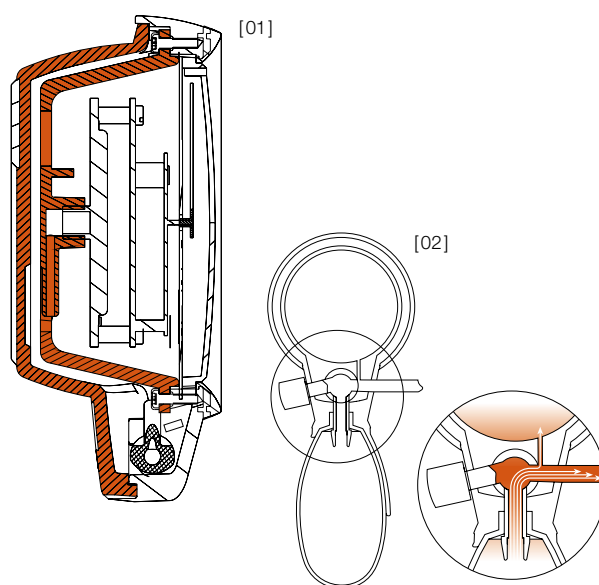

AIR PRESSURE DISTRIBUTION SYSTEM

- Single-tube construction with air pressure distribution system [02]



CLINICAL APPROVAL

HEINE GAMMA® G7

The model with aluminium alloy construction and precision valve



- ▶ **Shockproof**
- ▶ **Latex-free**
- ▶ **Unique pressure measuring system**
- ▶ **Pressure-proof manometer**
- ▶ **Micro filter.** Protects valve and manometer
- ▶ **Large diameter scale.** 56 mm diameter black-on-white scale is easy to read
- ▶ **High quality indicator.** Precisely calibrated, durable and shock resistant
- ▶ **Large insufflation bulb.** 50 % larger bulb volume for fast inflation of the cuff
- ▶ **Exclusive air-release valve.** Provides fine and precise air release rate adjustment, quick deflation, and extreme durability.
- ▶ **Upper housing made of aluminium alloy** with a shock absorbing ring protecting indicator and glass
- ▶ **Lower housing is moulded of shockproof plastic** with shock absorbing, anti-slip coating
- ▶ **Light alloy, adjustable spoon.** For both right- or left-handed operation
- ▶ **High quality, single-tube cuff with velcro fastener.** Washable
- ▶ **Ergonomic design.** G7 is perfectly balanced for comfort and ease of use
- ▶ **10 year guarantee***

* against defects in design, materials and workmanship for a period of 10 years from the date of purchase.

GAMMA G7

complete with adult cuff, in zipper pouch	M-000.09.232
complete with child cuff, in zipper pouch	M-000.09.233
practice kit with child, adult and small adult cuff in zipper pouch	M-000.09.554
complete with adult cuff without zipper pouch, box of 10	M-000.09.560

HEINE GAMMA® G5

The HEINE standard sphygmomanometer



Combines all design and quality features of the GAMMA G7, except:

- ▶ **Housing is moulded of shockproof plastic** with shock absorbing, anti-slip coating
- ▶ **Adjustable spoon (plastic).** For both right or left-handed operation
- ▶ **5 year guarantee***

* against defects in design, materials and workmanship for a period of 5 years from the date of purchase.

GAMMA G5

complete with adult cuff, in zipper pouch	M-000.09.230
complete with child cuff, in zipper pouch	M-000.09.231
practice kit with child, adult and small adult cuff in zipper pouch	M-000.09.555
complete with adult cuff without zipper pouch, box of 10	M-000.09.561

HEINE GAMMA GP

The HEINE standard sphygmomanometer with push-button valve



Combines all design and quality features of the GAMMA G5, but with a push-button valve:

► **Innovative push-button valve.** Fine and precise air-release rate adjustment, ultra-fast air-release, practically wear-free. Easy to operate by both left and right-handed users

► **5 year guarantee***

* against defects in design, materials and workmanship for a period of 5 years from the date of purchase.

GAMMA GP

complete with adult cuff, in zipper pouch	M-000.09.242
complete with child cuff, in zipper pouch	M-000.09.243
practice kit with child, adult and small adult cuff in zipper pouch	M-000.09.556
complete with adult cuff without zipper pouch, box of 10	M-000.09.562

HEINE GAMMA® XXL LF

The top-quality professional model

Since the GAMMA XXL LF is non-mercurial, it is both safe and environmentally-friendly.




[01]

- ▶ **Latex-free**
- ▶ **Durable polycarbonate construction.** Will hold up in the busiest hospital and practice environments
- ▶ **Oversized 130 mm diameter scale.** Can be pivoted for clear, easy reading from any angle
- ▶ **Long spiral inflation tube (3 m).** For increased working distance
- ▶ **High quality indicator.** Precisely-calibrated and durable
- ▶ **Precision air-release valve.** Fine and precise air release
- ▶ **Maintenance-free,** precise and durable manometer
- ▶ **Mercury-free.** Environmentally friendly
- ▶ **Wall, desk, rail and wheeled stand.** Adjustable height
- ▶ **10 year guarantee***

*against defects in design, materials and workmanship for a period of 10 years from the date of purchase.

GAMMA XXL LF

GAMMA XXL LF-W wall complete with adult cuff [01]	M-000.09.323
GAMMA XXL LF-S wheeled stand complete with adult cuff [02]	M-000.09.326
GAMMA XXL LF-R rail complete with adult cuff [03]	M-000.09.325
GAMMA XXL LF-T desk complete with adult cuff (non-spiral insufflation tube) [04]	M-000.09.322
GAMMA XXL LF-WB for wallboard, see page 159 [05]	M-000.09.327

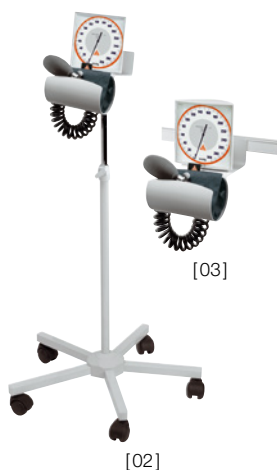

[02]

Accessories

cuff-holder for GAMMA XXL LF-T	M-000.09.141
desk stand for GAMMA XXL LF-T	M-000.09.142
adjustable wheeled stand for GAMMA XXL LF-S (height 80 – 125 cm)	M-000.09.145
universal mounting plate for GAMMA XXL LF-R, -S, -W	M-000.09.306
rail clamp for GAMMA XXL LF-R	M-000.09.307



[04]



[05]


Blood-pressure cuffs

Latex-free blood-pressure cuffs


Cuffs for GAMMA G7 | G5 | GP | XXL LF sphygmomanometers.




- ▶ **Latex-free**
- ▶ **Fit the most widely used sphygmomanometers**
- ▶ **High-quality velcro cuffs.** Physiologically tested
- ▶ **Double-sewn**
- ▶ **Washable**
- ▶ **Index mark for exact positioning**
- ▶ **Clear indication of size and control mark for arm circumference**
- ▶ **Extra thick, highly flexible tube.** For a long life
- ▶ **Proven reliability.** Absolutely airtight system

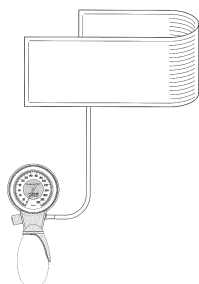

GAMMA G7 G5 GP	Single-tube cuffs (cpl.)		
	Size	Circumference cm	
infant [01]	6x28	8–13	M-000.09.700
child [02]	10x34	13–20	M-000.09.702
adult small [03]	14x53	20–29	M-000.09.704
adult [04]	14x58	29–41	M-000.09.619
adult large [05]	17x70	35–47	M-000.09.706
thigh [06]	20x86	41–59	M-000.09.615

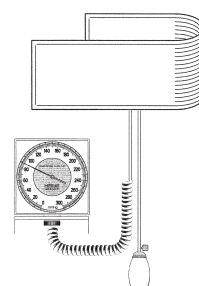

GAMMA XXL LF	Twin-tube cuffs (cpl.)		
	Size	Circumference cm	
infant [01]	6x28	8–13	M-000.09.701
child [02]	10x34	13–20	M-000.09.703
adult small [03]	14x53	20–29	M-000.09.705
adult [04]	14x58	29–41	M-000.09.629
adult large [05]	17x70	35–47	M-000.09.631
thigh [06]	20x86	41–59	M-000.09.625

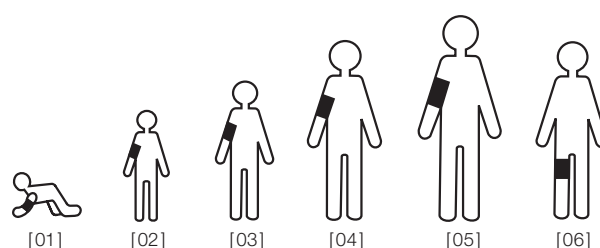

HEINE GAMMA® stethoscopes



- ▶ **Solid chest piece** and ideal weight avoid stray noise from manipulation and ensure excellent resonance. Hygienic matt chrome finish
- ▶ **Chest piece with non-chill rims for the comfort of the patients**
- ▶ **Ergonomic, chrome-plated binaurals.** Individually adjustable
- ▶ **High-quality Y-piece tubing** for perfect acoustics and sound reproduction. Length 56 cm
- ▶ **Latex-free materials**

GAMMA 3.1 pulse stethoscope For fine acoustic performance



- ▶ **Stethoscope with solid flat chest piece.** The chest piece can be easily inserted under the sphygmomanometer cuff
- ▶ **High-quality membrane,** 44 mm diameter
- ▶ **Single-tube construction** with excellent sound transmission
- ▶ **Includes ear olives** (large|soft) for sound insulation
- ▶ **Optional available:** ear olives (small|hard) (M-000.09.946)

GAMMA 3.1

M-000.09.941

GAMMA® 3.2 acoustic stethoscope

For excellent acoustic performance



Technical data as GAMMA 3.1, but with a double chest piece.

- ▶ **Combined double chest piece** with two different frequency ranges to choose from: A high-quality membrane on one side and a conical acoustic head on the other side
- ▶ **High-quality membrane**, 44 mm diameter for excellent acoustic performance and sound transmission
- ▶ **Cone diameter 30mm**

GAMMA 3.2	M-000.09.942
-----------	--------------

GAMMA 3.3 acoustic stethoscope

For the examination of children



Technical data as GAMMA 3.2

- ▶ **Reduced contact area** of 31 mm makes GAMMA 3.3 ideal for examining children

GAMMA 3.3	M-000.09.943
-----------	--------------